




GRAND CANYON  
FOUNDATION™

# GCF LIGHT

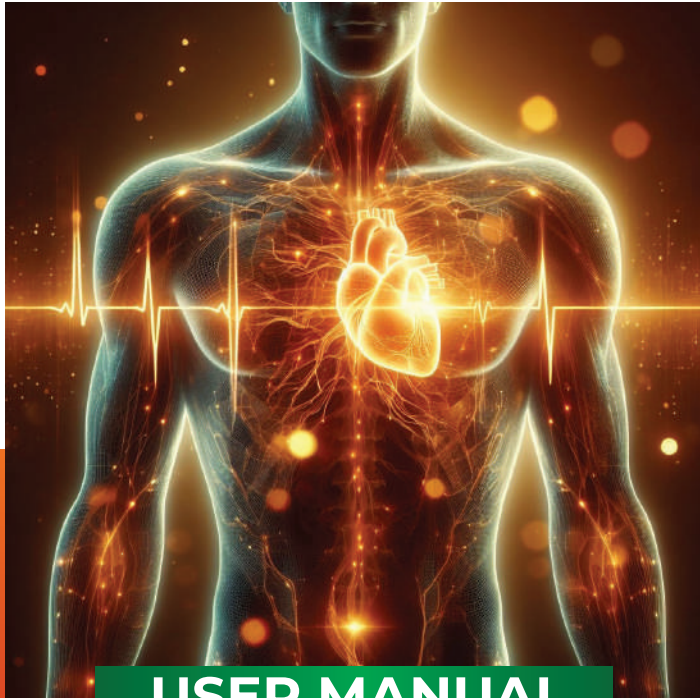

## USER MANUAL

Meet the future of personal wellness. The GCF Light is a powerhouse using precise wavelengths of light therapy together with Frequency Nutrient technology that is portable, safe and easy to use.



GRAND CANYON  
FOUNDATION™

# Table Of Contents

---

Included With This Package

---

Benefits

---

Safety First

---

How to Operate your GCF Light

---

Product Specification



GRAND CANYON  
FOUNDATION™

# Welcome **Message**



Congratulations on taking part in a new and exciting paradigm of self care! The GCF Light is the first biophotonic device that combines LED light waves with our revolutionary technology that transforms high quality nutrients into frequencies that the body can use.

This hand-held device uses a combination of 5 LEDs that emit light in the following therapeutic wavelengths: blue (460nm), red (630nm 660nm), and infrared (850nm & 940nm). In addition to the LEDs, 30 nutritional supplements are added in the form of frequencies.



GRAND CANYON  
FOUNDATION™

## Included With This **Package**

---

Zippered case  
The GCF Light  
Attachment  
USB-C charging cable  
Protective Eye Wear  
Cleaning Cloth





GRAND CANYON  
FOUNDATION™

# GCF Light **Benefits**

---

- Enhance circulation
- Reduce inflammation
- Support cellular repair
- Boost collagen production
- Address acne-causing bacteria
- Ease sore muscles
- Recover faster from workouts



GRAND CANYON  
FOUNDATION™

# Safety First

---

This device is not waterproof. Do not submerge the GCF light and its components in water.

The GCF light is automatically timed for 3 minutes, then beeps and shuts off. For continued use on a treatment area or multiple areas, allow at least 1 minute off to cool before starting again.

Do not stare directly at the light.

Do not treat close to eyes. Wear protective glasses while using the device.

Do not use if you have epilepsy.

Do not use if your skin is light sensitive or currently using medication or undergoing treatment that may cause skin to be sensitive to light.

Ask your doctor prior to using if you are pregnant or suspect you are pregnant as we don't have long term studies available at this time.

Consult with your doctor prior to using with your child.

Consult with your veterinarian prior to using with your pets.



GRAND CANYON  
FOUNDATION™

# How to Use **Your GCF Light**

---

## *Charging your GCF Light*

OPTION 1: Plug into Outlet

OPTION 2: Plug into an optional mobile power bank (Not included)

Note:

- Uses built in rechargeable 18650 battery
- Type C charging port
- Battery indicator shows battery status

## **ON/OFF Button**

Turn ON: press the power button and release. (Do not hold it down).

Turn OFF: device will beep and shut off automatically after 3 minutes use.

You may turn it off at any time by pressing and releasing the power button.



GRAND CANYON  
FOUNDATION™

# Light Wave Options



The GCF Light has 5 LEDs in the following light spectrum:  
BLUE (460 nm) RED (630 nm, 660 nm) INFRARED (850 nm, 940 nm)  
You can choose from 3 light wave options:

- (A) ALL 5 LEDs (red, blue infrared). Default Option: Press start button & release.
- (B) Red LEDs only: After pressing start button and release, press again & release
- (C) Blue LED only: After pressing start button and release, press 2 more times & release





GRAND CANYON  
FOUNDATION™

# How to Use **Modes**

---

The GCF Light offers both steady and pulsating light referred to as modes. A unique feature of the GCF light is the option to select the way light is delivered to the body. Our research shows that the pulsating mode provides enhanced absorption, deeper tissue repair, increased cellular activation, superior therapeutic results.

## **Steady Mode (default mode):**

A constant stream of light appears when you press the power button and release.

## **Pulsating Mode:**

To activate pulsating mode, press the power button, press again and hold.



GRAND CANYON  
FOUNDATION™


# How to Use Your GCF Light Attachment

To begin, screw the attachment onto the hand-held light.

Remove the red cap prior to use. We recommend replacing the red cap onto the attachment after each use.



The attachment is used for focusing on smaller areas of the body such as ears, nostrils, teeth and gums.

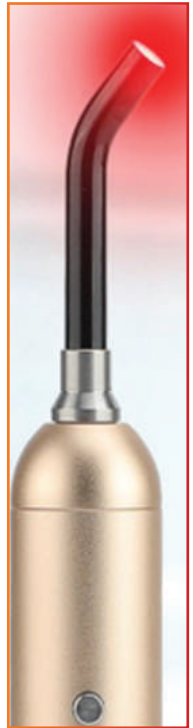



GRAND CANYON  
FOUNDATION™

# Product Specification

Rating Power	15W
Actual Output	5W
Product Size	145 mm * D30 mm
Product Weight	0.135 kg including one battery side
LED QTY	5 PCS LEDs
Wavelength	460nm 630nm 660nm 850nm 940nm
Material	Alumium for the torch body Glass for the attachment
Battery Type	18650 Lithium Rechargeable Battery
Casing Color	Champagne-Gold
Battery	Rechargeable Lithium Battery 18650
Irradiation	Over 120 mW / cmsq
EMF	0 UT at 0 Inch (SAFE)
Timer	3 Minutes
Beep Sounds	Yes
Lifespan	50,000 hours



# GRAND CANYON FOUNDATION™



Grand Canyon Foundation 10775 Double R Blvd. Ste 888  
Reno, NV 89521  
1 (800) 722-3870  
[support@grandcanyonfoundation.com](mailto:support@grandcanyonfoundation.com)  
[www.grandcanyonfoundation.com](http://www.grandcanyonfoundation.com)

---

These Claims Have Not Been Evaluated by The Food and Drug Administration.  
These Products are Not Intended To Diagnose, Treat, Cure or Prevent Any Disease.